P/N 800065 - 10

THD SLIDE ONE KIT



<u>Use</u>

THD Slide One Kit provides all instruments needed to perform the THD procedure; treatment of varying grades of haemorrhoidal disease.


General description

THD Slide is specifically designed for immediate use with the THD Revolution/THD Evolution doppler dearterialiser. THD Slide anoscope/proctoscope is specifically designed for improving the THD procedure while treating the higher grades of the haemorrhoidal disease. The specific shape of the anoscope/proctoscope, which includes a port for the THD Doppler Probe, allows ligation and prolapse reduction with a single tool and reduces procedure time. The THD Slide anoscope/proctoscope is made of clear, medical grade, plastic which creates, in tandem with the use of the fibre optic tip, ideal lightning of the operative field. The THD Slide anoscope/proctoscope is designed to maintain the correct tension of the anal cavity during the procedure and presents placement for the THD needle holder. The supplied doppler probe lodges into the anoscope/proctoscope and is thereby positioned with the anal mucosa for an accurate detection of the haemorrhoidal arteries. The THD Slide One box contains 10 individually packed disposable Doppler probes, which can be immediately be used without preventive cleaning and sterilisation. The box contains 10 adsorbable sutures packs (6 sutures each pack) and 10 Anal Spon haemostatic sponges.

Product Datasheet

P/N 800065 - 1 P/N 800065 - 5 P/N 800065 - 10

ΕN

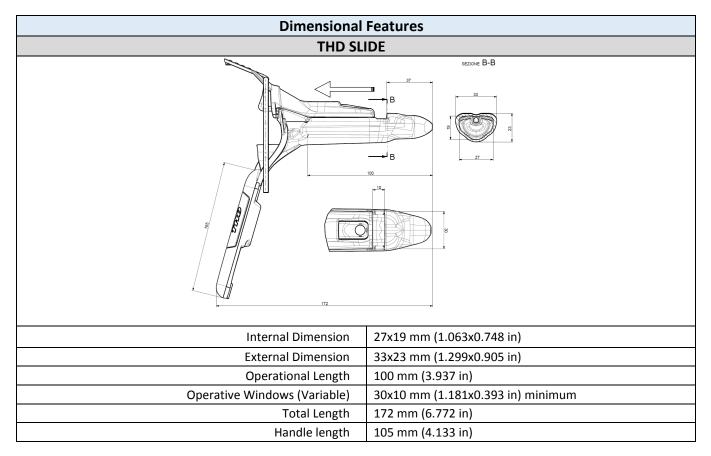

THD Slide One (Sterile)						
P/N		Q.ty	Description			
800065-10		10	THD Slide One Kit (sterile)			
include	710004	10	THD Slide (sterile) — with 1 proctoscope, 1 needle holder, 1 pusher			
include	800063	10	Disposable Doppler Probes (sterile)			
include	880000	10	n. 6 THD Suture (sterile)			
include	800050	10	Scrofalon Anal (sterile)			

THD Slide One (Sterile)					
P/N	Q.ty	Description			
800065-1	1	THD SLIDE ONE KIT			
800065-5	5	THD SLIDE ONE KIT			
800065-10	10	THD SLIDE ONE KIT			

Associated Equipment						
P/N	Qty	Description				
800133	1	THD REVOLUTION				
800000	1	THD EVOLUTION				

THD360025_EN_Rev01_20190326_THD Slide One.docx