



THD SURGY STERILE DISPOSABLE



plastic materials of medical grade and it specifically designed fibre optic bended tip

Use

The SURGY anoscope is a SINGLE-USE and STERILE medical device indicated for anorectal canal analysis and surgery.

General description

Its compact design allows an easy access to the anus, thus a good dilation and a gentle retraction of the mucosa, without any trauma to the sphincter. The wide lateral slit allows an optimal visualisation of the anal canal and enhanced the operative field.

The asymmetrical design stretches the anal canal laterally and it avoids the prolapse of the mucosa.

THD Surgy is made with high quality transparent

is designed to be connected to an Illuminator, through the


Features

- Made of medical-grade polycarbonate
- Complies with the biocompatibility requirements established by ISO 10993
- Can be used with THD curve fiber optics
- Atraumatic insertion profile
- Extremely compact
- Ready to use
- Lightweight (40g)
- Sterilised with Ethylen Oxide
- Class I CE device
- Individually packed
- Latex-free



Advantages:

- Sterile device
- Does not require preparation before use
- Eliminates the risk of cross-contamination
- Compact and easy to handle
- Mini-invasive and atraumatic

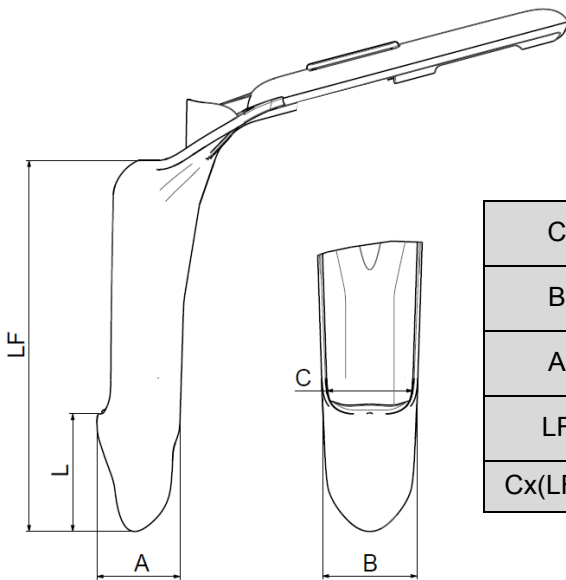

		Surgy mini	Surgy Maxi
C	Internal lumen	23 mm	28 mm
B	Outside dimension	26 mm	31 mm
A	Outside dimension	23 mm	27 mm
LF	Device length	122 mm	122 mm
Cx(LF-L)	Window	23x85 mm	28x85 mm

THD Surgy LightScope Sterile		
Code	Qty	Denomination
800004	10	THD Surgy Mini Sterile Disposable
800005	10	THD Surgy Maxi Sterile Disposable