




# THD LIGHTSCOPE RECTOSCOPE STERILE

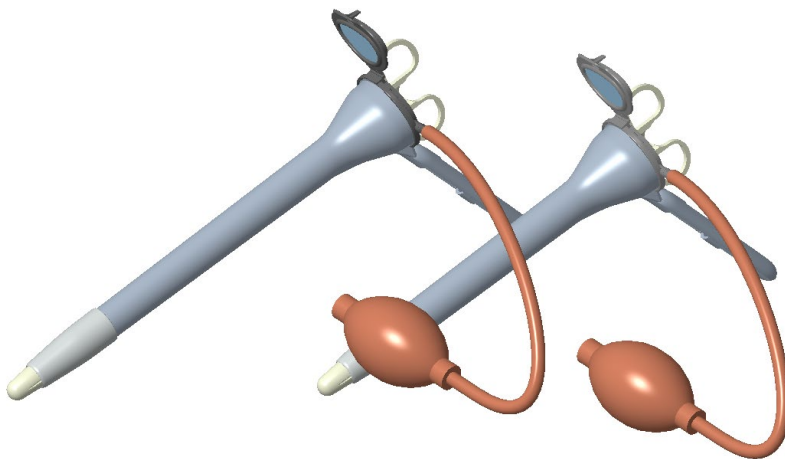

## Use

The rectoscope THD LightScope is a DISPOSABLE and STERILE medical device, indicated for analysis of ano-rectal canal and the distal region of the sigmoid.

## General Description

The rectoscope THD LightScope integrates in the cap an hand lens. The handle has the housing for use with

the light source THD LED2000 and THD REVOLUTION (light source function).

The THD Light-Scope rectoscope is ready-to-use and it is also available with the blowing pump included in the sales package. THD Light-Scope Recto is available with two lengths (200mm and 250mm).

## Features

- Made of a medical grade Polycarbonate
- Comply with the biocompatibility requirements of the ISO 10993 standard
- Two length available 200 mm and 250 mm (7.87 in and 9.84 in)
- Lens with a 4x magnification
- Available also with disposable insufflator
- Ready to use
- Light (55 g/0.12 pound)
- CE Class Is
- Individually packaged
- Latex Free
- Sterilized by Ethylene Oxide



# Technical Data

P/N. 800120  
P/N. 800121  
P/N. 800122  
P/N. 800123

EN

### Advantages

- It doesn't require a preparation before its use
- It delete the risk of cross-contamination
- Sterile Device

Dimensional Features		
THD LIGHTSCOPE RECTO STERILE		
	<b>Mini</b>	<b>Maxi</b>
Diameter (B)	23 mm (0.90 in)	23 mm (0.90 in)
Length (A)	200 mm (7.87 in)	250 mm (9.84 in)

THD LightScope Recto Sterile		
Code	Qty	Denomination
800120 - 10	10	THD LightScope Recto Maxi Sterile
800121 - 10	10	THD LightScope Recto Mini Sterile
800122 - 10	10	THD LightScope Recto Maxi with insufflator Sterile
800123 - 10	10	THD LightScope Recto Mini with insufflator Sterile