




THD PRESSPROBE

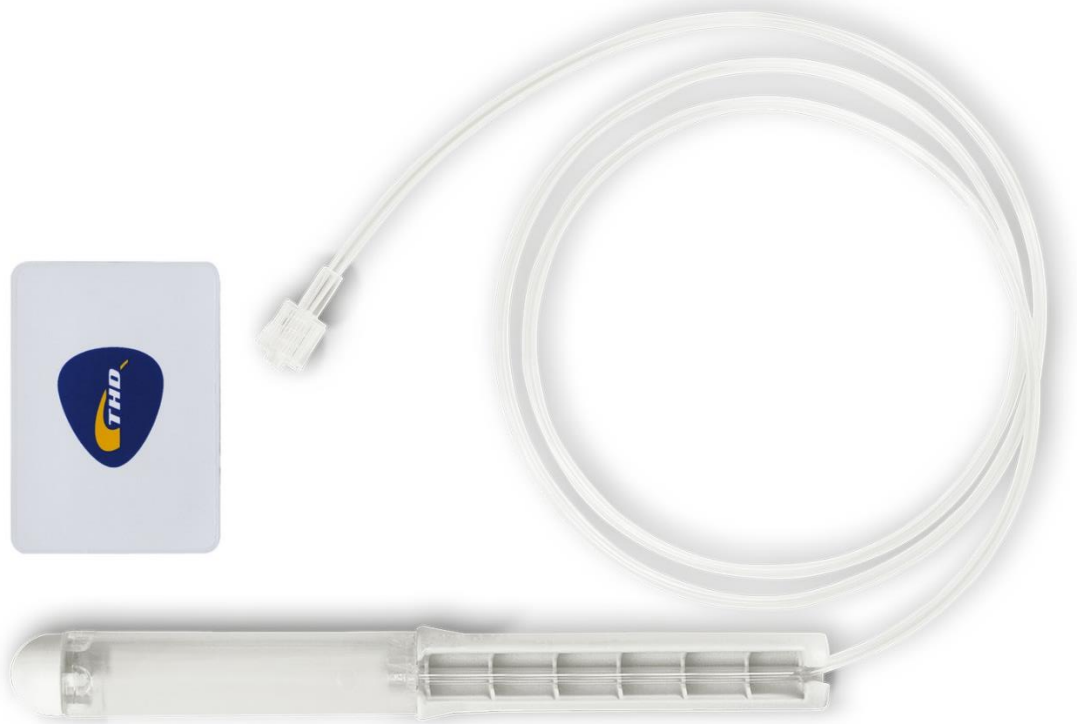

THD PressProbe is the innovative disposable active probe for the execution of the THD anal manometry exam.

The probe allows for the measurement of the real medium sphincter tone pressure in the different functional activations (basal, squeeze and strain tones).

THD PressProbe is meant to be used exclusively in combination with THD Anopress.

There is no need of calibrating the probe, because of the self-calibrating routine before the start of the procedure.

The probe has an ergonomic grip and a completely atraumatic profile, which ensures the maximum comfort for both doctor and patient during the examination and it has been designed to reproduce a normal digital examination.

Each pack contains a proprietary Probe Card, to enable the recording of the examination.



Probe Characteristics	
Probe Diameter	16 mm
Probe length	165 mm (including the handle)
Length of measuring space	50 mm
Length of tubing	1000 mm
Probe Weight	23 g
Connection to Anopress	Luer Lock male

PressProbe – PressKit*			
Code	Denomination	Q.ty	Description
800147	THD PRESSPROBE (BOX)	10	Pack of 10 PressProbes in a single box
800150	THD PRESSPROBE (ENV)	50	Pack of 50 PressProbes in a single plastic envelope
800148	THD PRESSPROBE KIT (TR BOX)	10	PressProbes and Anoscope N-ano truncated (box)
800151	THD PRESSPROBE KIT (TR ENV)	50	PressProbes and Anoscope N-ano truncated (envelope)
800149	THD PRESSPROBE KIT (FL BOX)	10	PressProbes and Anoscope N-ano bevelled (box)
800152	THD PRESSPROBE KIT (FL ENV)	50	PressProbes and Anoscope N-ano bevelled (envelope)

*see codes datasheet 800040, 800134, 800165 and 800166

CONNECTED DEVICES		
Code	Denomination	Description
800140	THD ANOPRESS	Anorectal manometry probe
800160	THD ANOPRESS	Anorectal manometry probe
800198	THD PROCTOSTATION	Integrated Proctologic Platform
800194	THD PROCTO MOBILE	Mobile Proctologic Platform